



DuoTEMP[®] 
dual cure temporary filling material

Double security keeps you smiling!

- Dual-cure
- High marginal integrity
- Easy removal
- One material with multiple indications

 **COLTENE**

DuoTEMP is a dual-cured, eugenol-free, white temporary filling material formulated from zinc oxide/zinc sulphate for optimal protection from the infiltration of bacteria. DuoTEMP is one of the first all-in-one temporary filling materials that can be cured using both light and saliva. This solution delivers greater flexibility, allowing the patient to immediately resume normal daily activities. DuoTEMP also contains a pleasant flavor that keeps the dentist and patient smiling.

► **Double security**

• **Dual-cure**

After only 40 s of light curing, the patient can immediately load the material. Afterwards, DuoTemp will continue to completely self-cure on its own with saliva.

• **High marginal seal due to slight expansion**

Reduces infiltration of bacteria

► **One material with multiple indications**

Reduces inventory and minimizes confusion

- Temporary fillings of class I and II cavities
- Temporary fillings for endodontic applications
- Temporary fillings for inlay / onlay preparations

► **Natural tooth-like feeling**

DuoTEMP's high surface hardness feels like a natural tooth to minimize irritation and maximize patient comfort.

► **Additional benefits**

- Contains zinc oxide and fluoride for additional protection
- Eugenol-free – no effect on the curing of composites
- Radiopaque

Order Information:

5830 DuoTEMP Intro Kit

2 syringes of 5 g, Intraoral Separator 5ml

5831 DuoTEMP Single Pack


1 syringe of 5 g

7090 DuoTEMP Eco Pack

5 syringes of 5 g

7235 Intraoral Separator, 5ml

► **Marginal Seal Comparison**

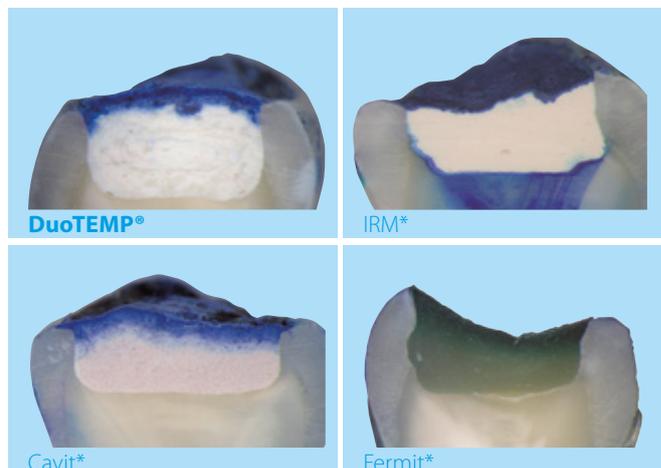

Cross-sectioned extracted human molars were filled with cured DuoTEMP and then immersed in methylene blue dye to determine the marginal sealing properties. DuoTEMP demonstrated excellent sealing properties.

Source: internal data

* Not a trademark of Coltène/Whaledent

► **Easy, time saving application**



One component without mixing. Good retention in the cavity. Easy removal.

© 2015 Coltène/Whaledent AG – www.coltene.com

Coltène/Whaledent AG

Feldwiesenstrasse 20
9450 Altstätten / Switzerland
Tel +41 71 757 5300
Fax +41 71 757 5301
info.ch@coltene.com

Coltène/Whaledent Inc.

235 Ascot Parkway
Cuyahoga Falls, Ohio 44223 / USA
Tel +1 330 916 8800
Fax +1 330 916 7077
info.us@coltene.com

Coltène/Whaledent Ltd.

The President Suite-A
Kendal House, Burgess Hill
Victoria Way
West Sussex, RH15 9NF/U.K.
Tel. +44 (0)1444 235486
Fax +44 (0)1444 870640
info.uk@coltene.com

